


FERNO TSL Extenders for Scoop®

Accessory to widen the patient surface on a Scoop® EXL for the transportation of larger patients

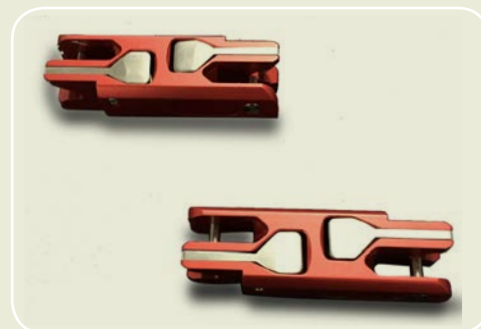

FOR MORE DETAILED
INFORMATION SEE OVER

Specifications

| Length | Width | Thickness | Load Limit | Weight |
|--------|-------|-----------|-------------------|-------------|
| 130mm | 35mm | 25mm | 227kg/36st/500lbs | 245g (each) |

FERNO – is a world leading innovator and manufacturer of patient transportation solutions.

Our product range is born of a closer understanding of our customers' needs and the practical demands of their work. Ferno's industry-leading **TSL Extenders** typifies our passion for top quality, highly functional design... patient transfer products that don't just do the job, they do it **BETTER**.

call +44 (0) 1274 851 999
email sales@ferno.co.uk
www ferno.co.uk

 **FERNO**®

FERNO TSL Extenders for Scoop®

Accessory for widening the surface area of Scoop® EXL for the transportation of physically larger patients

With a safe working load of 227 kg (36 stone) and adjustable length, Ferno's Scoop® EXL is capable of carrying a range of patients, from young to old, small to large. For the physically larger-size patients, for example bariatric patients and sports casualties such as rugby players, larger body masses and shapes sometimes cause issue in terms of physical support and manouverability on the Scoop®.

When attached to the head-end, foot-end or at both ends, the Ferno TSL Extender accessory allows the Scoop® to be extended in shape and width to accommodate people of differing body shapes and sizes.

Combined with the impressive 227 kg SWL, this allows the Scoop® to be used for the majority of patients that require transferring and transporting.

- **Easy to apply** – TSL extenders simply clip in between each half of the Scoop®
- **Flexibility of use** – can be used at head-end, foot-end or at both ends
- **Accommodate patients of all sizes** – use one or 2 extenders according to the patient's size and shape
- **Augment the use of existing equipment** - can be used on existing or newly purchased Scoops®



Part Nos:

25-0030-001 TSL Extender (single) for Scoop® EXL

FERNO - TSL extenders FPC001 02.16

call +61 (7) 3881 4999
email info@ferno.com.au
www ferno.com.au



Ferno Australia, 11 Johnstone Road, Brendale, Queensland 4500