



VERSATILITY MEETS URGENCY.

The Ferno XTRIC8 Rescue Sled is a universal extrication solution engineered for speed, reliability, and adaptability across diverse rescue environments.

Whether you're navigating lowland terrain, scaling mountainous regions, entering confined spaces, or executing high-angle extrications, XTRIC8 delivers unmatched performance when every second counts.


The XTRIC8 is a compact, durable rescue sled designed for rapid deployment. Its intelligent design reduces the time and effort required to secure casualties, while retaining traditional functionality and tactical evacuation capabilities.


Key Features

- » Supplied as a complete, ready-to-use unit with all essential hardware and lifting bridles included.
- » Pre-rigged for immediate deployment out of the bag.
- » An integrated internal patient harness ensures optimal casualty security during vertical lifts and high-angle extractions.
- » Colour coded, 4-point lifting bridle is pre-set for rapid and efficient implementation
- » Reflective external webbing for enhanced visibility in low light conditions
- » Straps and buckles are glove-friendly and designed for efficient attachment and release
- » Robust carry bag with reflective webbing, four large pockets and adjustable backpack harness
- » XTRIC8 is intuitive to use and reduces length of user training



XTRIC8 Load Rating

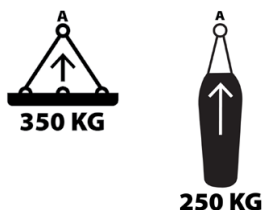

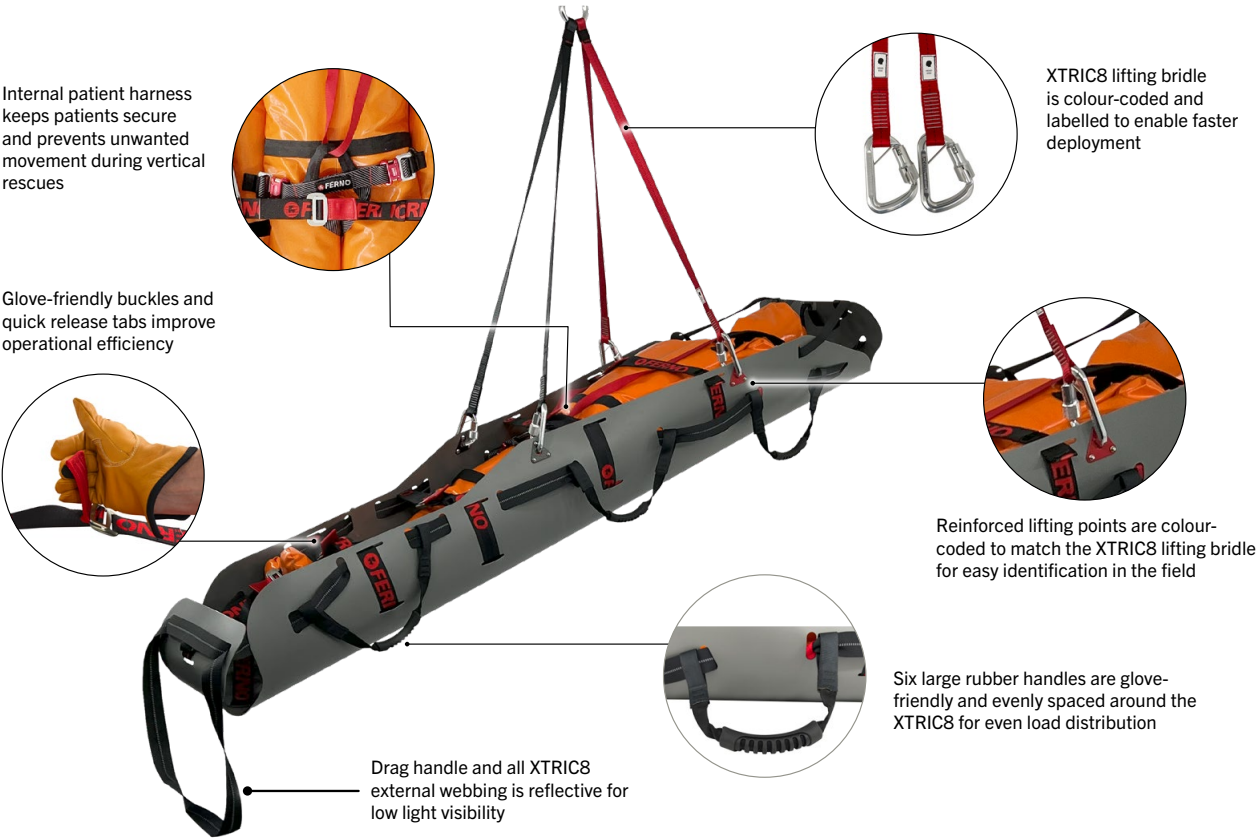

SPECIFICATIONS

XTRIC8 RESCUE SLED

Product Code	FWE XTRIC8
Inclusions	XTRIC8 Rescue Sled pre-rigged with internal patient harness assembly, external restraints, external straps and handles
	1 x Maillon Rapide for vertical lifting
	1 x Colour coded 4-Point Lifting Bridle
	1 x Storage Bag with backpack straps
Length	2540 mm
Width - Widest Point	914 mm
Weight - XTRIC8 only (rigged)	7 kg
Weight - XTRIC8 (rigged) and Bag	8.5 kg
Load Limit - XTRIC8	350 kg - HORIZONTAL / 250 kg - VERTICAL
XTRIC8 LIFTING BRIDLE	
Weight - Lifting Bridle	0.9 kg
Load Limit - Lifting Bridle	500 kg
Bridle Karabiners	VYO N-2853G - Aluminium Alloy (SAE 7075) Triple Action Karabiner
Bridle Karabiner Dimensions	124 mm x 80 mm / Gate Opening: 22 mm
Bridle Karabiner Ratings	Minimum Tensile Strength: 23 kN / Minimum Gate Strength: 16 kN
Bridle Karabiner Standards	ANSI Z359.12-2019; CSA Z259.12-16; EN 362:2004

- Ferno and manufacturers of components reserve the right to change specifications without notice.
 - General specifications are rounded to the nearest whole number. Metric conversions are calculated before rounding the English measurements.
- For more information contact Ferno Customer Service or your local distributor.